



E1000 Dx Digital Pathology Solution

Designed for your team's unique workflow needs
from digitisation to sign out



Epredia's comprehensive solution complements your existing workflow, increasing team satisfaction and efficiency

E1000 Dx Scanner

Maximise productivity, minimise workflow disruptions and capture the best image the first time.

Maximise throughput

- **Dual Slide Processing:** Eliminates downtime between sides, as the instrument can immediately start scanning the next slide
- **Continuous Loading:** Load slides on demand directly from your coverslipper and prioritise slides based on your workflow needs
- **Safety Container:** Avoids scanning disruption by segmenting slides that need additional attention, while continuing the digitisation run

Automated Quality Control

- Gain confidence that your workflow will not be interrupted due to poor quality images
- Automated focus quality algorithm rescans slides with advanced settings, if out-of-focus areas are detected
- Proprietary sample-detection algorithm, produces high quality focus and image resolution

True Walkaway Scanning

- Remotely monitor your scanners throughput data and workload statistics
- Email alerts throughout the scanning process for run completion and any errors
- Overnight runs of up to 1000 slides, add an additional shift to your laboratory

Integrated IMS and Viewer Software

Integrated image management software and viewer, that allows seamless movement and review of images for caseload allocation and review.

Easy to Navigate

- Case priority feature, allows you to enhance your case review based on urgency

Facilitates Collaboration

- Telepathology capable allowing collaboration between internal and external colleagues
- Integration with LIS system allows you to ensure all critical data is visible in the IMS and viewer for each case and patient

Efficient Review Process


- Nine-slide view: Open and link up to nine slides to quickly navigate multiple slides from the case
- Browser compatibility gives pathologists freedom to review cases from any location on their iPhone or iPad



Continuous
Flexible
Intuitive

A Scalable and Secure Solution

Providing easy access for pathologists on site or in remote locations

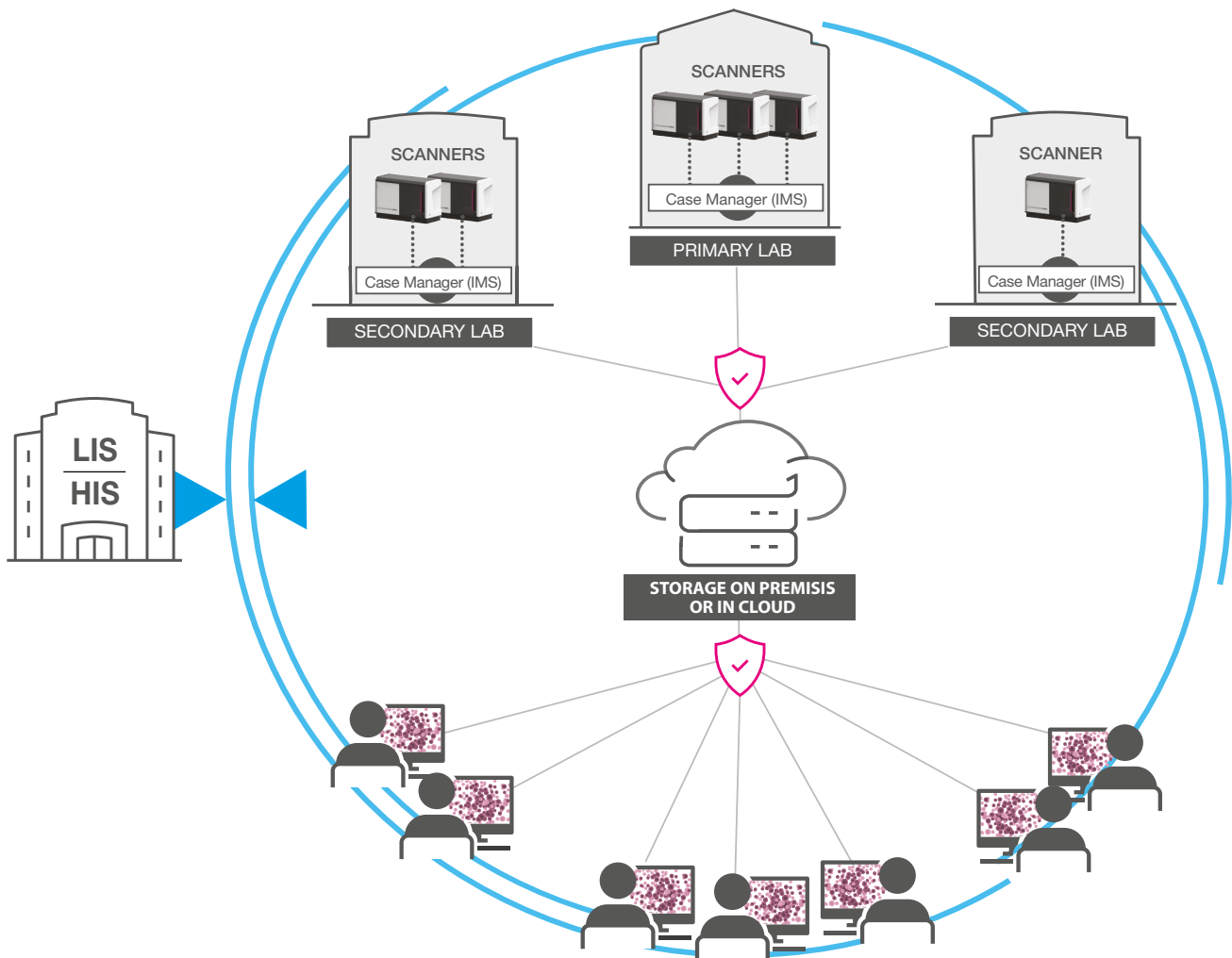

The Epredia Experience.

At the heart of successful pathology labs.

Powered by people passionate about your success, The Epredia Experience gives you increased control of your total workflow with a holistic approach that combines trusted solutions and extraordinary service.

The Epredia Experience is different.
It makes a difference for you and your lab.



Product Specifications

Capacity: 1000 slides

Loading: Continuous and Batch

Throughput: 72 slides per hour

Scan Speed: 28 seconds

Dimensions: (WxDxH):
154X100X91CM

Weight: 270 kg

Throughput based on 15X15
sample at 40x resolution scan

Find out more at epredia.com

Tudor Road | Manor Park | Runcorn | WA7 1TA

© 2023 Epredia Group. All rights reserved. All trademarks are the property of Epredia Group and its subsidiaries unless otherwise specified. M2400027-UK_A4_R032024

epredia
A Member of PHC Group